




# SpecView®/4848 Controller

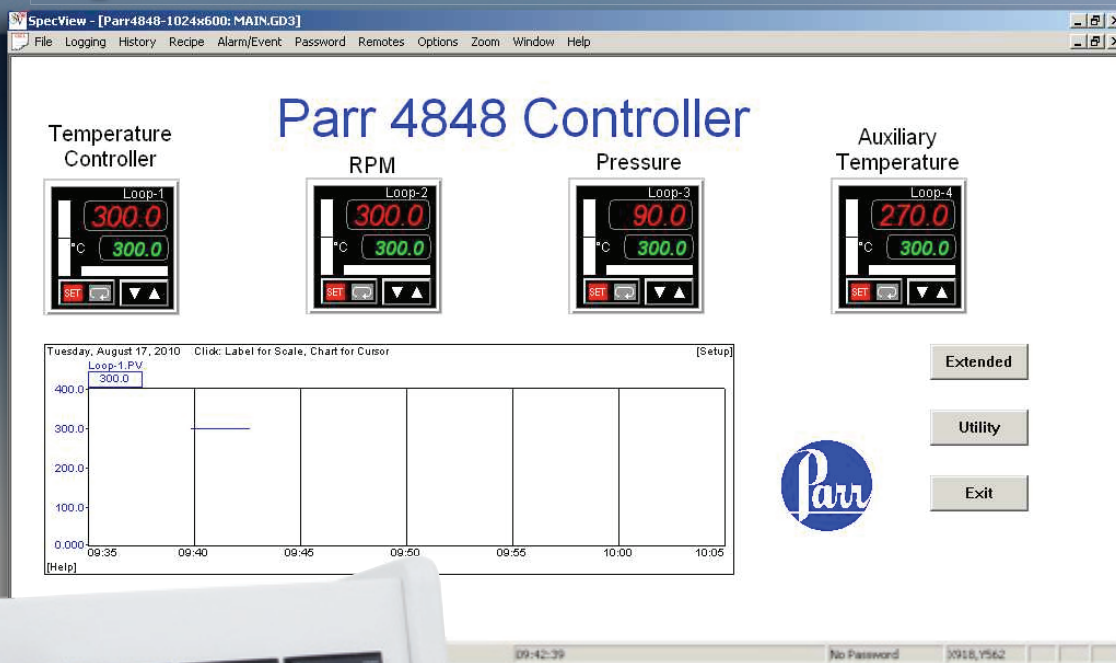

Sample screenshot from a PC running SpecView Software with data from a Parr® 4848 Process Controller.

## Features and Benefits

In comparison to the standard ParrComm software, this new SpecView build offers the following enhancements:

- Improved real-time plotting including multiple charts each with an unlimited number of pens
- Enhanced data logging — including the ability to log any controller parameter, not just the process variable
- Alarm and event logging
- Superior alarm annunciation
- The ability to easily customize the user interface to a given application
- The user interface can be localized, for example to support a variety of languages.

## Ordering Guide

Order No.	Description
A3504HC	Parr 4848/SpecView Package

Each SpecView Package Includes:

2181E	Parr 4848/SpecView License Dongle
578M	SpecView Software Installation CD
577M	SpecView Operator Instruction Manual

Contact Parr Instrument Company or any of its Dealers/Agents for more information, including prices and availability.



For nearly a decade, Parr Instrument Company has used SpecView as the principal means to provide a flexible, customizable user interface for the 4870 Series Process Controllers. The unified Modbus communications architecture of the 4848 Controller system allows us to now bring the power and convenience of SpecView to the 4848 Series Controller. A special arrangement with the developers of SpecView allows Parr to offer a unique version of SpecView with support for 4848 Series instruments incorporating up to eight individual loop or limit controllers. This provides support for up to two fully populated 4848 Controllers or one 4848B Controller.